C3 Male Continence Device

With a unique design that fits your anatomy


A Discreet and Effective way to Prevent or Decrease Urine Leakage



Comfort with C3

The contoured, cradle-like design and elastic strap of the C3 minimizes excessive pressure by distributing the radial force evenly around the penis. Carefully chosen materials also enhance your comfort... soft foam wrapped in supple plastic.

Available in two sizes, the C3 offers a customized fit for even greater comfort.

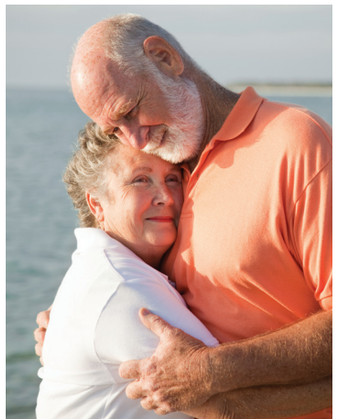

Control with C3

When the C3 is placed on the penis, a urethral occlusion pad presses on the urethra. While this localized pressure prevents the involuntary loss of urine, the flow of blood is unimpeded to assure your safety.

Confidence with C3

So compact, the C3 is easy to conceal. The elastic strap holds the device in position, allowing you to go about your daily routine without concern about the C3 shifting.

Proper hygiene never has to be a concern either. The disposable C3 is reasonably priced, so you can replace the device at regular intervals. The C3 is a more economical choice than absorbent products.



U. S. Patent No. 5,184,629 and related global patents. C3 is a registered trademark of Personal Medical Corp.

© 1998-2012 Personal Medical Corp. All rights reserved.

290046-001 Rev N

personalmed
advancing pelvic health®

8672 154th Ave. NE
Redmond, WA 98052
Tel 866.839.9260 • +1 425.242.6373
Fax +1 425.497.1045
www.personalmed4men.com

C3® Male Continence Device



A Solution That Offers Comfort, Control and Confidence

personalmed
advancing pelvic health®

C3 Male Continence Device

Patient Order Form

Patient, please complete form and mail with payment to:

Personal Med
8672 154th Ave. NE
Redmond, WA 98052

Name

Address

City

State Zip Code

Telephone

____ Bags of 6 Units @ $60.00 = ____

+ Shipping & Handling = $9.00
(Continental US only)

+ Sales Tax (WA only) = ____

Total Amount of Order = ____
(Payment Required with Order)

The C3 Male Continence Device is sent to you in discreet packaging.

Questions? Call 866-839-9260

Method of Payment (Please do not send cash):

❑ Check payable to Personal Med. Allow 14 days to process.

❑ Money Order/Cashier's Check

❑ Visa ❑ Master Card ❑ American Express

Card Number Expiration Date

Signature Security Code

Prescribing Physician Information form must accompany this form.

C3 Male Continence Device

Prescribing Physician Information

Must be completed by physician.

Personal Med
8672 154th Ave. NE
Redmond, WA 98052

866-839-9260
425-242-6373
Fax: 425-497-1045

Patient Name

Physician Name

Physician Address

City State

Zip Code Telephone

Specialty UPIN#

Principle Diagnosis

ICD-9 Code Number

Size Around: ❑ Small less than 2½ inches ❑ Regular 2½-4 inches ❑ Large 4-5½ inches

❑ Permanent urinary incontinence #Refills ____
(necessary for Medicare filing)

Contributing Diagnoses

	Description	ICD-9 Code No.
Primary Cause		
Secondary Cause		

Medicare HCPCS Billing Code: A4356
Description: External Urethral Clamp or Compression Device

Questions? Call 866-839-9260

I prescribe the C3 Male Continence Device as non-surgical incontinence management for this patient.


Physician Signature Date

The C3 Male Continence Device must be prescribed by a physician. Please ask your doctor to complete and sign the Prescribing Physician Information form.



Easy to Use

Here's how it works: Unfasten the strap on a flat device and place your penis through the opening. The foam side should face in and the strap should be on top.

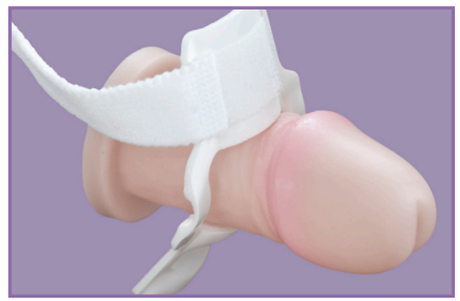

Fold the device towards your body to close it around your penis. (The hinges where the device sticks out should be AWAY from the body and the section with the bump should be underneath your penis.) Position the device an inch or so away from your body before tightening the strap.



Choose the Size that's Right for You

The C3 device comes in three sizes. Follow these instructions to determine the size that's right for you:

1. Use a flexible tape measure.
2. Measure only a non-erect penis.
   For best results, conduct sizing when your body is at room temperature. For example, don't conduct the sizing immediately after a shower.
3. Wrap the tape measure *around* your penis, one to two inches from the body. (Do not measure the length.)
4. When the tape measure is snug around your penis, read the number in inches. Do not pull tightly.
5. Use the number of inches around to determine your size before ordering:

   Small = less than 2½ inches (6 cm) around
   Regular = 2½-4 inches (6-10 cm) around
   Large = 4-5½ inches (10-14 cm) around

personalmed
advancing pelvic health

8672 154th Ave. NE
Redmond, WA 98052
866.839.9260
+1 425.242.6373
Fax +1 425.497.1045
www.personalmed4men.com